

**GETINGE WD15 CLARO
THE CLEAR CHOICE IN
WASHER-DISINFECTORS**





INSIGHT IS WHAT MAKES THE DIFFERENCE

Every time we develop a new product, we base our design work on our deep insight into our customers' everyday reality. It's a successful approach, and it has made us the world leader in sterilization and disinfection.

When we designed the new Getinge WD15 Claro our focus was clear: how can we create the optimal under-counter washer-disinfector for a dental practice? How can we ensure the best possible efficiency, ergonomics and ease-of-use?

The result is a high-capacity, fast-working and nice-looking machine called Claro.

The nickname comes from one of its most distinguishing feature, which actually represents another kind of "insight" – that you have a clear view of what happens throughout the decontamination process, thanks to the glass window and the chamber light.

ALL THE CAPACITY YOU NEED IN 600 MM

Claro is perfectly suited for dental practices where a compact under-counter machine is often the ideal choice to save countertop space. Of course, it combines very well with Getinge tabletop sterilizers and sealers and there is also a wide range of practical accessories.



Unique chamber design

Thanks to smart positioning of water/air connections, the chamber has a large capacity. The carts offer an almost 100% usable surface.

Clever self-closing system

Thanks to the clever self-closing mechanism for the spray arm docking system, the upper cart, along with its spray arm, can be removed to make space for larger goods, keeping the same level of cleaning efficiency.

By a simple click the spray arm is completely opened at both ends simplifying the maintenance of the Claro.





VIP – View-In-Process glass window

For extra safety you have a constant clear view of the process thanks to the glass window and the lighting. It's easy to see by just a glance that the spray wings rotate as they should and that the goods are loaded correctly.

Integrated hot air drying

After rinsing, washing and disinfection phases, there is an extra built-in drying phase with HEPA 13 filter that ensures that the goods are quickly and efficiently dried.

Pump pressure monitoring

For extra safety, the circulation pump has a pressure sensor that stops the process and generates an error code if the water pressure is not correct.

Easy to operate

The PACS 350 control system is very user-friendly and program settings are easily changed on the operator panel.

Documentation via printer or Netcom

Printer or Netcom with USB connection for detailed process reports is available as an option. Claro can also be connected to T-DOC, Getinge's comprehensive instrument management system, via Netcom.

Getinge Online

Getinge Online gives you insight into the status and performance of your Getinge washer-disinfectors. You can get secure access from any computer or smart phone anywhere in the world. With access to comprehensive and detailed information, process printouts, archive of errors and warnings (including instructions for how to resolve them), automatically updated overviews of all your Getinge equipment and much more. And with Getinge Online you can even setup if you or somebody from your team should get a text/SMS or e-mail when your washer-disinfector has completed a process or generated a message.

Clear display

The machine is self-diagnostic. Program, phase, chamber temperature, error codes and remaining cycle time are clearly shown in the display. It's easy to follow the process and it simplifies servicing.

Condenser cooling

Thanks to a smart solution with condenser cooling there is no need to connect to a ventilation system, which makes it less humid in the room and also simplifies installation.

Well-tested

Getinge's range of detergents has been developed, tested and validated on Getinge washer-disinfectors to consistently meet EN ISO 15883 standards for cleaning efficiency of goods commonly found in CSSD's, clinics and practices to protect users and patients.



ALL THE ACCESSORIES YOU MAY NEED

Getinge offers a well-adapted range of accessories and there's a solution whatever kind of goods you need to reprocess.



LOWER BASIC CART

With support for larger loads.
Capacity: 2 SPRI II (340x250x70 mm) or 2 DIN 1/1 (480x250x50).
Free loading height: with upper wash cart 320 mm.
Without upper wash cart 490mm. Loading area 520x496mm.
Art.no. 503014600



UPPER BASIC CART

Capacity: 2 SPRI II (340x250x70 mm) or 2 DIN 1/1 (480x250x50).
Free loading height: 110 mm.
Art.no. 503013400




2 LEVEL WASH CART

Capacity: 4 SPRI II (340x250x70 mm) or 4 DIN 1/1 (480x250x50 mm)
Free loading height: 1 level to spray wing 120 mm, spray arm 90 mm
2 level 80 mm
Loading area: 510x470 mm
Art.no. 503393600
Not possible to use this cart together with any injection bar



DENTAL INJECTION BAR

With 10 injectors for hand pieces. Five injectors are for hand pieces with an approximate outer diameter of 16 mm, five for 20 mm. A fine mesh filter is included. (Hand pieces shown in pictures are not included)
Art.no. 503043200

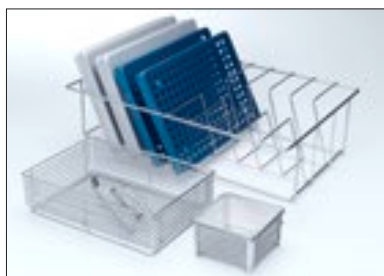

DENTAL KIT

With 1 tray rack for 9 dental trays, 1 SPRI III instrument tray and 1 fine mesh basket with lid for small parts.
(Trays on pictures not included)
Art. no: 503707200

AVAILABLE IN TWO COLORS

The Getinge WD15 Claro is available in blue/stainless steel or in white powder coated. Both variants can be delivered with a cabinet for a more ergonomic work height.



Getinge WD15 Claro

Design	Blue or White
External Dimensions	596 x 850 x 625 mm (W x H x D)
Internal Dimensions	538 x 520 x 534 mm (W x H x D)
Weight	92 kg
Chamber volume	180 L
Effective chamber volume	150 L
Chamber/Panels Material	AISI 316L/AISI 304L/galvanized
Communication Ports	RS232 and RS485
Electrical	Single-or 3-phase 50/60 Hz 200V, 208V, 220-230V, 240V or 400V
Heating Element	Water 6.3 kW dryer 2.8 kW
Water Supply	Cold and Hot
Dosage Pumps	2
Detergent Dosage Control	Yes (standard)
Circulation Pump	380 L/min
Noise	<60 dB (A)
Standards	EN ISO 15883, MDD, EN IEC 61010-1, EN 61326-3-2, DVGW
Optional	Cabinet adjustable from 468-520 mm 2 extra dosing pumps with flow control Water softener External printer NetCom bluebox T-DOC (instrument management system)



Getinge Infection Control is the world leading provider of solutions for sterile processing in the healthcare sector. We aim to ensure the highest quality and safety at the lowest total cost. We offer complete solutions for a seamless work-flow, reducing the risk of contamination while helping healthcare to increase efficiency.

Getinge Infection Control AB

P O Box 69
SE-305 05 Getinge
SWEDEN
Phone: +46 (0) 10 335 00 00
E-mail: info@getinge.com

GETINGE
GETINGE GROUP

getinge.com

GETINGE GROUP is a leading global provider of products and systems that contribute to quality enhancement and cost efficiency within healthcare and life sciences. We operate under the three brands of ArjoHuntleigh, GETINGE and MAQUET. **ArjoHuntleigh** focuses on patient mobility and wound management solutions. **GETINGE** provides solutions for infection control within healthcare and contamination prevention within life sciences. **MAQUET** specializes in solutions, therapies and products for surgical interventions, interventional cardiology and intensive care.